

ComLab Front Panel Valve 1/8" Inlet/Outlet for Quick Connect Fittings - Dry Media- AIR/ CO2/ HE/ AR/ N2

Enware laboratory front panel valves are designed for fitting to the front panels or mullions of Fume Cupboards or similar equipment. The durable powder coat finish gives significant protection against corrosion, discolouration and surface damage. Slimline handles are ergonomically designed for precise feel and fingertip control. Handles are colour coded for easy service identification and use complimentary buttons. Designed to be connected to incoming service and remote outlet nozzle.



KEY FEATURES

- Suited for dry media up to technical gas level 4.0
- Chemically and ultrasonically cleaned
- Can be fixed from the front with screws (not supplied) or from the rear with backnut
- Teflon needle valve with stainless steel seat
- Durable white epoxy powder coat finish with coloured metal handles
- 1/8" BSP FI Inlet/Outlet for use with Quick Connect or Compression Fittings
- Blue handle and indicator button identification to DIN EN13792
- Indicator buttons for compatible media included
- Compatible media: Compressed Air, Carbon Dioxide, Helium, Argon, Nitrogen

PRODUCT CODE



LCN150NRDRY	ComLab Front Panel Valve 1/8" Inlet/Outlet for Quick Connect Fittings - Dry Media- AIR/ CO2/ HE/ AR/ N2
-------------	---

Due to ongoing Research and Development, specifications may change without notice.
Component specifications may change on some export models.
Refer to warranty statement for warranty details - www.enware.com.au/warranty.

ComLab Front Panel Valve 1/8" Inlet/Outlet for Quick Connect Fittings - Dry Media- AIR/ CO2/ HE/ AR/ N2


TECHNICAL DATA

Finish	Powder Coating
Valve Type	Needle Valve
Mounting Type	Front Panel Valve
Maximum Recommended Pressure	700kPa
Maximum Testing Pressure	1000kPa
Inlet Connection	1/8" BSP FI
Outlet Connection	1/8" BSP FI
Connection Notes	Requires quick connect fittings LCN173 or Compression Fittings
Gas Options	Argon, Carbon Dioxide, Compressed Air, Helium, Nitrogen

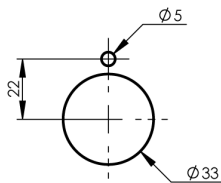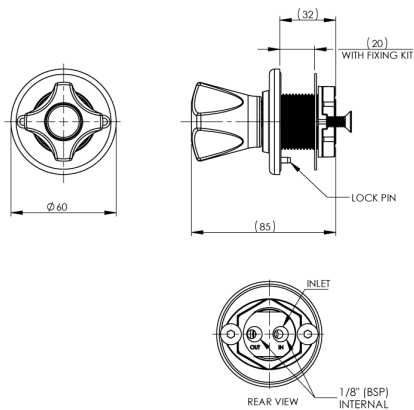

INSTALLATION DETAIL
(MOUNTING SURFACE CUT OUT)

PLEASE CHECK LOCKING PIN
POSITIONS BEFORE INSTALLATION
AS ORIENTATION MAY VARY

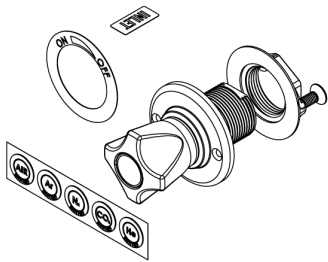

REQUIRED OPTIONS

LCN173	Spare Parts: Quick Connect Straight - 1/8" BSP x 8mm
--------	--

ComLab Front Panel Valve 1/8" Inlet/Outlet for Quick Connect Fittings - Dry Media- AIR/ CO2/ HE/ AR/ N2

SPARE PARTS

LCN184	Spare Parts: ComLab - Stuffing Box Assembly - Needle Valve - LCN Products
LCN191-BLUE	ComLab: Handle (Blue)
LCN193	ComLab: Indicator Mount
LCN181	ComLab: Button Replacement - Dry Media - LCN Products
LCN173	Spare Parts: Quick Connect Straight - 1/8"BSP x 8mm